

SCP SCIENCE presents the

UNIVERSAL SCRUBBER SYSTEM

Remove and Capture
more than 98% of
all fumes generated !!



All Acid Resistant
Tubing !!

"Take your digestion system out of the fume hood"

The Universal Scrubber System (USS) uses a 3 vessel array in combination with a robust vacuum system to scrub out up to 98% of the acid fumes generated by an acid digestion. The USS is perfect for either TKN or trace metal digestion fumes. Adapters are available to fit the USS onto either type of digestion system. Free up fumehood space and stop using "water-jet" vacuum with the **SCP SCIENCE** Universal Scrubber System.


Components

- 1- Three 4-liter polypropylene bottles complete with;
 - a) A diffusion system for increased fume neutralization
 - b) PVC adaptors and caps
- 2- Vacuum Pump System
- 3- Speciality cap to drain and rinse bottles
- 4- Connector for *DigiVAC* evacuation system
- 5- Acid resistant tubing and PVC connections

Ordering Information & Specifications

UNIVERSAL SCRUBBER SYSTEM

Catalog Number	010-520-060
Vacuum	42 CFM draw (at 0" H ₂ O) 33 CFM (at 10" H ₂ O)
Sound rating	55 db-c
Power	115V/230V (User Selection)
Size - Bottle rack	8.5" x 25"
- Pump	12" x 12"

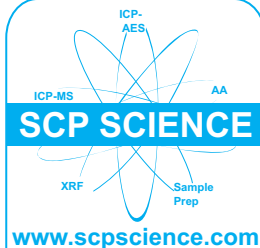

Sample
Prep

ICP-AES & MS

AA

XRF

Calibration
Standards

Certified
Reference
Materials

Use the Universal
Scrubber System
on your *DigiPREP HT*
High Temperature
Digestion System !!



The *DigiPREP HT* is
a new System for
Kjeldahl digestion

- Versatile heating to 450°C
- Acid resistant graphite block
- Available in two sizes:
40 tube (100ml) or
20 tube (250ml)
sample capacity

Call for more information
or to request a quotation!!

USA

Tel.: (800) 361-6820
Fax: (800) 253-5549



Canada / International

Tel.: (800) 361-6820 / (514) 457-0701
Fax: (800) 253-5549 / (514) 457-4499



Europe

Tel.: +33 (0)1 69 18 71 17
Fax: +33 (0)1 60 92 05 67