


CERTIFICATE OF ANALYSIS FOR

Copper-Cobalt Ore (Democratic Republic of the Congo) CERTIFIED REFERENCE MATERIAL OREAS 550

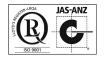

Document: COA-1475-OREAS 550-R1

TABLE OF CONTENTS

INTRODUCTION	8
SOURCE MATERIAL	8
PERFORMANCE GATES	8
COMMINUTION AND HOMOGENISATION PROCEDURES	9
PHYSICAL PROPERTIES	9
ANALYTICAL PROGRAM	9
STATISTICAL ANALYSIS	10
Homogeneity Evaluation	16
PARTICIPATING LABORATORIES	16
PREPARER AND SUPPLIER	17
METROLOGICAL TRACEABILITY	22
COMMUTABILITY	22
INTENDED USE	22
STABILITY AND STORAGE INSTRUCTIONS	23
INSTRUCTIONS FOR CORRECT USE	23
HANDLING INSTRUCTIONS	25
LEGAL NOTICE	25
DOCUMENT HISTORY	25
QMS ACCREDITATION	25
CERTIFYING OFFICER	25
REFERENCES	25
APPENDIX	26
LIST OF TABLES	
Table 1. Certified Values and Performance Gates for OREAS 550	3
Table 2. Indicative Values for OREAS 550.	7
Table 3. Physical properties of OREAS 550.	9
Table 4. 95% Confidence & Tolerance Limits for OREAS 550	11
LIST OF FIGURES	
Figure 1. Co by 4-acid digestion in OREAS 550	18
Figure 2. Co by Peroxide Fusion in OREAS 550	19
Figure 3. Cu by 4-acid digestion in OREAS 550	20
Figure 4. Cu by Peroxide Fusion in OREAS 550	21

Table 1. Certified Values and Performance Gates for OREAS 550.

							Relative Standard Deviations 5% windo				
Constituent	Certified		Absolute	Standard	Deviations	3	Relative	Standard D	eviations	5% w	indow
Constituent	Value	1SD	2SD Low	2SD High	3SD Low	3SD High	1RSD	2RSD	3RSD	Low	High
Borate Fusion	XRF										
Al ₂ O ₃ , wt.%	5.39	0.059	5.27	5.51	5.21	5.57	1.09%	2.19%	3.28%	5.12	5.66
BaO, ppm	386	33	319	453	285	486	8.68%	17.36%	26.04%	366	405
CaO, wt.%	9.70	0.115	9.47	9.93	9.35	10.04	1.18%	2.37%	3.55%	9.21	10.18
Co, wt.%	0.145	0.004	0.138	0.152	0.135	0.156	2.44%	4.88%	7.32%	0.138	0.152
Cu, wt.%	0.278	0.011	0.257	0.299	0.247	0.310	3.78%	7.57%	11.35%	0.264	0.292
Fe ₂ O ₃ , wt.%	1.89	0.022	1.85	1.94	1.83	1.96	1.14%	2.28%	3.42%	1.80	1.99
K ₂ O, wt.%	0.709	0.012	0.686	0.732	0.674	0.743	1.62%	3.25%	4.87%	0.674	0.744
MgO, wt.%	10.72	0.086	10.55	10.89	10.47	10.98	0.80%	1.60%	2.39%	10.19	11.26
MnO, wt.%	0.160	0.001	0.157	0.163	0.156	0.164	0.87%	1.73%	2.60%	0.152	0.168
P ₂ O ₅ , wt.%	0.081	0.005	0.071	0.092	0.066	0.097	6.43%	12.86%	19.29%	0.077	0.085
SiO ₂ , wt.%	52.79	0.319	52.16	53.43	51.84	53.75	0.60%	1.21%	1.81%	50.15	55.43
SO ₃ , wt.%	0.565	0.018	0.529	0.601	0.511	0.619	3.19%	6.39%	9.58%	0.537	0.593
SrO, ppm	< 100	IND	IND	IND	IND	IND	IND	IND	IND	IND	IND
TiO ₂ , wt.%	0.314	0.008	0.298	0.329	0.290	0.337	2.53%	5.06%	7.59%	0.298	0.329
Thermogravin	netry										
LOI ¹⁰⁰⁰ , wt.%	16.91	0.142	16.62	17.19	16.48	17.33	0.84%	1.68%	2.52%	16.06	17.75
Infrared Comb	oustion										
C, wt.%	4.11	0.105	3.90	4.32	3.80	4.43	2.55%	5.09%	7.64%	3.91	4.32
S, wt.%	0.220	0.012	0.195	0.245	0.183	0.257	5.61%	11.22%	16.84%	0.209	0.231
Sulphuric Acid	d 5% Leach										
Co, wt.%	0.101	0.017	0.067	0.135	0.050	0.152	16.74%	33.48%	50.23%	0.096	0.107
Cu, wt.%	0.197	0.012	0.172	0.221	0.160	0.233	6.22%	12.43%	18.65%	0.187	0.207
Sulphuric Acid	d 10% Leac	h									
Co, wt.%	0.097	0.015	0.068	0.127	0.053	0.142	15.19%	30.37%	45.56%	0.092	0.102
Cu, wt.%	0.201	0.024	0.154	0.249	0.130	0.272	11.73%	23.46%	35.20%	0.191	0.211
Peroxide Fusi	on ICP										
Ag, ppm	< 1	IND	IND	IND	IND	IND	IND	IND	IND	IND	IND
Al, wt.%	2.87	0.066	2.74	3.00	2.67	3.07	2.31%	4.63%	6.94%	2.73	3.01
B, ppm	229	14	202	256	188	269	5.92%	11.84%	17.76%	217	240
Ba, ppm	349	12	325	372	314	384	3.33%	6.66%	9.99%	331	366
Bi, ppm	1.45	0.29	0.87	2.03	0.57	2.33	20.16%	40.32%	60.48%	1.38	1.52
Ca, wt.%	7.03	0.136	6.75	7.30	6.62	7.43	1.94%	3.87%	5.81%	6.67	7.38
Ce, ppm	39.3	2.31	34.7	44.0	32.4	46.3	5.87%	11.74%	17.62%	37.4	41.3
Co, wt.%	0.150	0.005	0.140	0.159	0.135	0.164	3.18%	6.37%	9.55%	0.142	0.157
Cs, ppm	0.60	0.059	0.49	0.72	0.43	0.78	9.76%	19.53%	29.29%	0.57	0.63
Cu, wt.%	0.279	0.009	0.261	0.298	0.251	0.307	3.31%	6.63%	9.94%	0.265	0.293
Dy, ppm	3.57	0.300	2.97	4.17	2.67	4.47	8.39%	16.77%	25.16%	3.39	3.75
Er, ppm	1.93	0.21	1.51	2.34	1.31	2.55	10.72%	21.43%	32.15%	1.83	2.02
Eu, ppm	0.92	0.062	0.79	1.04	0.73	1.10	6.77%	13.55%	20.32%	0.87	0.96
Fe, wt.%	1.33	0.044	1.24	1.41	1.19	1.46	3.33%	6.65%	9.98%	1.26	1.39
Ga, ppm	7.80	0.748	6.30	9.30	5.56	10.05	9.59%	19.18%	28.77%	7.41	8.19
Gd, ppm	4.14	0.364	3.41	4.86	3.05	5.23	8.79%	17.59%	26.38%	3.93	4.34
Ho, ppm	0.66	0.08	0.49	0.82	0.41	0.90	12.37%	24.73%	37.10%	0.62	0.69
SI unit equival	antai nam	oorto nor	million = I	ma/ka = 11		001 14 0/	- 1000 pp	h norto no	r billion		·

SI unit equivalents: ppm, parts per million \equiv mg/kg \equiv μ g/g \equiv 0.0001 wt.% \equiv 1000 ppb, parts per billion.

Note 1: intervals may appear asymmetric due to rounding.



		Absolute Standard Deviations Relative Standard Deviation								5% w	indow
Constituent	Certified		1	1	1	1	Relative		I I	3 70 W	ıııdow
	Value	1SD	2SD Low	2SD High	3SD Low	3SD High	1RSD	2RSD	3RSD	Low	High
Peroxide Fus	ion ICP con	tinued									
K, wt.%	0.628	0.029	0.569	0.687	0.539	0.716	4.70%	9.40%	14.10%	0.596	0.659
La, ppm	22.8	2.03	18.8	26.9	16.7	28.9	8.88%	17.77%	26.65%	21.7	24.0
Li, ppm	105	6	92	118	86	124	6.12%	12.23%	18.35%	100	110
Lu, ppm	0.26	0.022	0.21	0.30	0.19	0.32	8.78%	17.55%	26.33%	0.24	0.27
Mg, wt.%	6.44	0.098	6.24	6.63	6.14	6.73	1.52%	3.03%	4.55%	6.12	6.76
Mn, wt.%	0.123	0.005	0.113	0.133	0.108	0.137	4.00%	8.00%	12.00%	0.117	0.129
Nb, ppm	9.37	1.63	6.10	12.63	4.47	14.26	17.43%	34.86%	52.30%	8.90	9.83
Nd, ppm	19.4	2.2	14.9	23.8	12.6	26.1	11.58%	23.16%	34.73%	18.4	20.3
P, wt.%	0.038	0.005	0.027	0.049	0.022	0.054	14.18%	28.36%	42.54%	0.036	0.040
Pb, ppm	16.8	4.3	8.1	25.5	3.8	29.9	25.80%	51.61%	77.41%	16.0	17.7
Pr, ppm	5.13	0.443	4.24	6.01	3.80	6.46	8.65%	17.30%	25.95%	4.87	5.38
Rb, ppm	22.9	1.68	19.5	26.3	17.9	28.0	7.35%	14.69%	22.04%	21.8	24.1
S, wt.%	0.221	0.018	0.184	0.257	0.166	0.276	8.26%	16.52%	24.79%	0.210	0.232
Sb, ppm	< 1	IND	IND	IND	IND	IND	IND	IND	IND	IND	IND
Si, wt.%	25.30	0.624	24.05	26.55	23.43	27.17	2.47%	4.93%	7.40%	24.04	26.57
Sm, ppm	4.37	0.396	3.58	5.16	3.18	5.56	9.08%	18.15%	27.23%	4.15	4.59
Sr, ppm	59	3.7	51	66	48	70	6.34%	12.69%	19.03%	56	62
Tb, ppm	0.61	0.038	0.53	0.68	0.49	0.72	6.21%	12.41%	18.62%	0.58	0.64
Th, ppm	5.37	0.392	4.58	6.15	4.19	6.54	7.30%	14.59%	21.89%	5.10	5.64
Ti, wt.%	0.188	0.009	0.169	0.206	0.160	0.215	4.88%	9.76%	14.64%	0.178	0.197
Tm, ppm	0.27	0.05	0.17	0.37	0.12	0.42	18.82%	37.64%	56.46%	0.25	0.28
U, ppm	23.9	2.9	18.2	29.6	15.3	32.5	11.97%	23.95%	35.92%	22.7	25.1
V, ppm	56	7	42	70	35	77	12.48%	24.96%	37.44%	53	59
Y, ppm	17.8	1.41	15.0	20.6	13.6	22.0	7.92%	15.84%	23.77%	16.9	18.7
Yb, ppm	1.70	0.170	1.36	2.04	1.19	2.21	9.99%	19.98%	29.97%	1.62	1.79
Zn, ppm	51	12	27	76	15	88	23.57%	47.13%	70.70%	49	54
*4-Acid Diges	tion										
Ag, ppm	0.127	0.018	0.092	0.162	0.075	0.180	13.78%	27.55%	41.33%	0.121	0.134
AI, wt.%	2.84	0.066	2.70	2.97	2.64	3.03	2.33%	4.66%	6.99%	2.69	2.98
As, ppm	5.82	0.59	4.63	7.00	4.04	7.59	10.16%	20.32%	30.48%	5.53	6.11
Ba, ppm	352	11	330	374	319	385	3.13%	6.26%	9.39%	334	369
Be, ppm	1.31	0.096	1.12	1.50	1.02	1.60	7.36%	14.71%	22.07%	1.25	1.38
Bi, ppm	1.40	0.086	1.23	1.58	1.14	1.66	6.15%	12.30%	18.45%	1.33	1.47
Ca, wt.%	6.91	0.217	6.47	7.34	6.26	7.56	3.15%	6.29%	9.44%	6.56	7.25
Cd, ppm	0.13	0.02	0.10	0.17	0.08	0.18	13.42%	26.84%	40.26%	0.13	0.14
Ce, ppm	36.6	3.08	30.4	42.7	27.3	45.8	8.41%	16.83%	25.24%	34.7	38.4
Co, wt.%	0.148	0.008	0.131	0.165	0.123	0.173	5.69%	11.38%	17.08%	0.140	0.155
Cr, ppm	23.1	6.6	9.8	36.3	3.2	42.9	28.64%	57.27%	85.91%	21.9	24.2
Cs, ppm	0.58	0.033	0.52	0.65	0.48	0.68	5.64%	11.28%	16.92%	0.55	0.61
Cu, wt.%	0.277	0.010	0.256	0.298	0.246	0.308	3.77%	7.54%	11.31%	0.263	0.291
Dy, ppm	3.30	0.107	3.09	3.52	2.98	3.62	3.22%	6.45%	9.67%	3.14	3.47
	lonto: nnm	1									

SI unit equivalents: ppm, parts per million \equiv mg/kg \equiv μ g/g \equiv 0.0001 wt.% \equiv 1000 ppb, parts per billion.



^{*}Four acid digestion quantitatively dissolves nearly all minerals in the majority of geological samples however, some refractory minerals may only be partially digested.

Note 1: intervals may appear asymmetric due to rounding.

		Absolute Standard Deviations Relative Standard Deviations								5% window	
Constituent	Certified		1	1	Г	I .	Relative	Otanuaiu D	- Eviations	3 70 W	iiidow
	Value	1SD	2SD Low	2SD High	3SD Low	3SD High	1RSD	2RSD	3RSD	Low	High
*4-Acid Diges	tion continu	ıed									
Er, ppm	1.76	0.128	1.50	2.01	1.37	2.14	7.27%	14.54%	21.81%	1.67	1.84
Eu, ppm	0.84	0.030	0.78	0.90	0.75	0.93	3.55%	7.11%	10.66%	0.79	0.88
Fe, wt.%	1.28	0.058	1.16	1.40	1.10	1.45	4.57%	9.14%	13.71%	1.21	1.34
Ga, ppm	7.65	0.285	7.08	8.22	6.79	8.50	3.73%	7.46%	11.19%	7.26	8.03
Gd, ppm	4.02	0.249	3.52	4.52	3.28	4.77	6.20%	12.40%	18.60%	3.82	4.23
Hf, ppm	2.19	0.180	1.83	2.55	1.65	2.73	8.23%	16.46%	24.69%	2.08	2.30
Ho, ppm	0.62	0.047	0.53	0.71	0.48	0.76	7.49%	14.98%	22.47%	0.59	0.65
In, ppm	0.18	0.013	0.15	0.20	0.14	0.21	7.09%	14.17%	21.26%	0.17	0.19
K, wt.%	0.614	0.034	0.545	0.682	0.511	0.716	5.57%	11.14%	16.71%	0.583	0.644
La, ppm	19.2	1.84	15.5	22.9	13.7	24.7	9.56%	19.12%	28.68%	18.3	20.2
Li, ppm	105	7	92	118	85	125	6.29%	12.58%	18.87%	100	110
Lu, ppm	0.23	0.016	0.20	0.27	0.19	0.28	6.65%	13.30%	19.95%	0.22	0.25
Mg, wt.%	6.43	0.229	5.97	6.88	5.74	7.11	3.56%	7.12%	10.67%	6.10	6.75
Mn, wt.%	0.119	0.006	0.108	0.131	0.102	0.137	4.98%	9.95%	14.93%	0.113	0.125
Mo, ppm	2.25	0.135	1.98	2.52	1.85	2.66	6.00%	12.00%	18.00%	2.14	2.37
Na, wt.%	0.040	0.009	0.022	0.058	0.013	0.067	22.61%	45.22%	67.83%	0.038	0.042
Nb, ppm	5.98	1.32	3.34	8.61	2.02	9.93	22.05%	44.10%	66.15%	5.68	6.28
Nd, ppm	18.5	0.81	16.8	20.1	16.0	20.9	4.38%	8.76%	13.13%	17.5	19.4
Ni, ppm	16.7	1.20	14.3	19.1	13.1	20.2	7.17%	14.35%	21.52%	15.8	17.5
P, wt.%	0.036	0.003	0.031	0.042	0.028	0.045	7.67%	15.34%	23.00%	0.034	0.038
Pb, ppm	13.7	0.90	11.9	15.5	11.0	16.4	6.57%	13.14%	19.72%	13.0	14.4
Pr, ppm	4.61	0.272	4.07	5.15	3.79	5.43	5.90%	11.80%	17.70%	4.38	4.84
Rb, ppm	22.7	1.45	19.8	25.6	18.4	27.1	6.39%	12.78%	19.17%	21.6	23.9
Re, ppm	0.008	0.001	0.005	0.010	0.004	0.012	17.51%	35.01%	52.52%	0.007	0.008
S, wt.%	0.233	0.014	0.205	0.261	0.191	0.275	6.00%	11.99%	17.99%	0.222	0.245
Sb, ppm	0.32	0.028	0.26	0.37	0.23	0.40	8.89%	17.77%	26.66%	0.30	0.33
Sc, ppm	6.27	0.310	5.65	6.89	5.34	7.20	4.94%	9.89%	14.83%	5.95	6.58
Sm, ppm	4.09	0.335	3.42	4.76	3.08	5.09	8.20%	16.40%	24.60%	3.88	4.29
Sn, ppm	1.18	0.13	0.92	1.45	0.79	1.58	11.17%	22.33%	33.50%	1.13	1.24
Sr, ppm	56	2.9	51	62	48	65	5.15%	10.29%	15.44%	54	59
Tb, ppm	0.62	0.035	0.55	0.69	0.52	0.73	5.64%	11.27%	16.91%	0.59	0.66
Th, ppm	5.15	0.461	4.23	6.07	3.77	6.53	8.95%	17.90%	26.84%	4.89	5.41
TI, ppm	0.094	0.008	0.077	0.111	0.069	0.119	9.01%	18.03%	27.04%	0.089	0.099
Tm, ppm	0.25	0.018	0.22	0.29	0.20	0.31	7.15%	14.31%	21.46%	0.24	0.26
U, ppm	24.7	2.9	19.0	30.4	16.1	33.2	11.56%	23.12%	34.69%	23.4	25.9
V, ppm	53	4.5	44	62	40	66	8.46%	16.93%	25.39%	50	56
W, ppm	0.78	0.16	0.47	1.10	0.32	1.25	19.87%	39.74%	59.61%	0.75	0.82
Y, ppm	16.5	1.10	14.3	18.7	13.2	19.8	6.69%	13.37%	20.06%	15.7	17.3
Yb, ppm	1.62	0.087	1.44	1.79	1.36	1.88	5.38%	10.75%	16.13%	1.54	1.70
Zn, ppm	41.9	2.46	37.0	46.8	34.5	49.3	5.88%	11.75%	17.63%	39.8	44.0
Zr, ppm	76	5.0	66	86	61	91	6.53%	13.07%	19.60%	72	80
	lonto: nnm				l	l	l .			L	L

SI unit equivalents: ppm, parts per million \equiv mg/kg \equiv μ g/g \equiv 0.0001 wt.% \equiv 1000 ppb, parts per billion.



^{*}Four acid digestion quantitatively dissolves nearly all minerals in the majority of geological samples however, some refractory minerals may only be partially digested.

Note 1: intervals may appear asymmetric due to rounding.

	Cartifical		Absolute	Standard		s	Relative	Standard D	eviations	5% w	indow
Constituent	Certified Value	1SD	2SD Low	2SD High	3SD Low	3SD High	1RSD	2RSD	3RSD	Low	High
Aqua Regia D	igestion										
Ag, ppm	0.086	0.010	0.065	0.107	0.055	0.117	12.04%	24.08%	36.12%	0.082	0.090
Al, wt.%	1.55	0.103	1.34	1.75	1.24	1.86	6.66%	13.32%	19.97%	1.47	1.63
As, ppm	4.93	0.304	4.32	5.54	4.02	5.84	6.17%	12.33%	18.50%	4.68	5.18
Ba, ppm	297	8	280	313	271	322	2.82%	5.64%	8.46%	282	311
Be, ppm	0.92	0.065	0.79	1.05	0.73	1.12	7.05%	14.11%	21.16%	0.88	0.97
Bi, ppm	1.07	0.043	0.98	1.15	0.94	1.20	4.03%	8.06%	12.09%	1.01	1.12
Ca, wt.%	6.70	0.190	6.32	7.08	6.13	7.27	2.84%	5.68%	8.52%	6.37	7.04
Cd, ppm	0.13	0.011	0.11	0.15	0.10	0.16	8.11%	16.22%	24.33%	0.12	0.14
Ce, ppm	15.8	1.33	13.1	18.5	11.8	19.8	8.44%	16.87%	25.31%	15.0	16.6
Co, wt.%	0.139	0.007	0.126	0.153	0.119	0.160	4.93%	9.86%	14.79%	0.132	0.146
Cr, ppm	21.5	1.30	18.9	24.1	17.6	25.4	6.04%	12.08%	18.12%	20.4	22.6
Cs, ppm	0.13	0.03	0.06	0.20	0.03	0.23	26.75%	53.50%	80.25%	0.12	0.14
Cu, wt.%	0.274	0.008	0.258	0.291	0.250	0.299	3.02%	6.05%	9.07%	0.261	0.288
Dy, ppm	1.92	0.076	1.77	2.08	1.70	2.15	3.94%	7.87%	11.81%	1.83	2.02
Er, ppm	0.85	0.040	0.77	0.93	0.73	0.97	4.74%	9.47%	14.21%	0.81	0.89
Eu, ppm	0.56	0.037	0.48	0.63	0.45	0.67	6.62%	13.24%	19.87%	0.53	0.58
Fe, wt.%	1.07	0.056	0.96	1.18	0.90	1.24	5.25%	10.49%	15.74%	1.02	1.13
Ga, ppm	4.45	0.296	3.86	5.04	3.57	5.34	6.64%	13.28%	19.92%	4.23	4.68
Gd, ppm	2.70	0.28	2.15	3.25	1.88	3.53	10.19%	20.38%	30.58%	2.57	2.84
Hf, ppm	0.31	0.06	0.19	0.42	0.13	0.48	18.98%	37.96%	56.94%	0.29	0.32
Hg, ppm	0.022	0.008	0.007	0.038	0.000	0.045	34.24%	68.48%	102.72	0.021	0.023
Ho, ppm	0.34	0.030	0.28	0.40	0.25	0.43	8.90%	17.79%	26.69%	0.32	0.36
In, ppm	0.16	0.007	0.14	0.17	0.13	0.18	4.49%	8.97%	13.46%	0.15	0.16
K, wt.%	0.126	0.015	0.097	0.156	0.082	0.171	11.75%	23.51%	35.26%	0.120	0.133
La, ppm	5.31	0.82	3.67	6.95	2.85	7.77	15.44%	30.88%	46.31%	5.05	5.58
Li, ppm	86	7.6	71	101	63	109	8.86%	17.73%	26.59%	82	90
Lu, ppm	0.099	0.008	0.083	0.116	0.075	0.124	8.31%	16.62%	24.93%	0.094	0.104
Mg, wt.%	6.08	0.201	5.68	6.48	5.48	6.69	3.30%	6.60%	9.90%	5.78	6.39
Mn, wt.%	0.118	0.008	0.103	0.133	0.095	0.141	6.40%	12.81%	19.21%	0.112	0.124
Mo, ppm	2.12	0.115	1.89	2.35	1.77	2.47	5.43%	10.86%	16.29%	2.01	2.23
Nd, ppm	8.32	0.523	7.27	9.36	6.75	9.88	6.28%	12.57%	18.85%	7.90	8.73
Ni, ppm	15.4	1.13	13.1	17.6	12.0	18.8	7.35%	14.71%	22.06%	14.6	16.1
P, wt.%	0.033	0.002	0.029	0.036	0.028	0.038	5.11%	10.22%	15.33%	0.031	0.034
Pb, ppm	11.9	1.10	9.7	14.1	8.6	15.2	9.21%	18.41%	27.62%	11.3	12.5
Pr, ppm	1.60	0.27	1.07	2.14	0.80	2.40	16.65%	33.30%	49.95%	1.52	1.68
Rb, ppm	4.16	0.50	3.15	5.17	2.64	5.67	12.14%	24.28%	36.42%	3.95	4.37
Re, ppm	0.007	0.001	0.005	0.009	0.004	0.010	12.45%	24.89%	37.34%	0.007	0.008
S, wt.%	0.224	0.019	0.185	0.262	0.166	0.281	8.56%	17.12%	25.68%	0.212	0.235
Sb, ppm	0.19	0.03	0.13	0.25	0.10	0.28	15.95%	31.90%	47.85%	0.18	0.20
Sc, ppm	4.44	0.327	3.79	5.10	3.46	5.42	7.36%	14.72%	22.08%	4.22	4.66
Se, ppm	1.11	0.25	0.61	1.62	0.36	1.87	22.70%	45.40%	68.09%	1.06	1.17
Sn, ppm	0.35	0.07	0.21	0.49	0.15	0.56	19.52%	39.04%	58.56%	0.33	0.37
Sr, ppm	44.3	2.94	38.4	50.1	35.4	53.1	6.64%	13.29%	19.93%	42.1	46.5

SI unit equivalents: ppm, parts per million \equiv mg/kg \equiv μ g/g \equiv 0.0001 wt.% \equiv 1000 ppb, parts per billion.

Note 1: intervals may appear asymmetric due to rounding.



Constituent	Certified		Absolute	Standard	Deviations	6	Relative	Standard D	eviations	5% w	indow
Constituent	Value	1SD	2SD Low	2SD High	3SD Low	3SD High	1RSD	2RSD	3RSD	Low	High
Aqua Regia Digestion continued											
Tb, ppm	0.39	0.012	0.36	0.41	0.35	0.42	3.14%	6.28%	9.42%	0.37	0.41
Te, ppm	0.035	0.010	0.015	0.056	0.004	0.066	29.26%	58.53%	87.79%	0.034	0.037
Th, ppm	2.23	0.186	1.85	2.60	1.67	2.78	8.38%	16.75%	25.13%	2.11	2.34
Ti, wt.%	< 0.01	IND	IND	IND	IND	IND	IND	IND	IND	IND	IND
TI, ppm	0.050	0.005	0.041	0.059	0.036	0.063	9.19%	18.38%	27.57%	0.047	0.052
U, ppm	22.3	1.74	18.8	25.8	17.1	27.5	7.82%	15.65%	23.47%	21.2	23.4
V, ppm	33.1	1.73	29.6	36.5	27.9	38.3	5.24%	10.47%	15.71%	31.4	34.7
W, ppm	0.28	0.06	0.17	0.39	0.11	0.44	19.75%	39.49%	59.24%	0.27	0.29
Y, ppm	8.81	0.460	7.89	9.73	7.43	10.19	5.22%	10.44%	15.66%	8.37	9.25
Yb, ppm	0.71	0.043	0.63	0.80	0.58	0.84	5.99%	11.98%	17.96%	0.68	0.75
Zn, ppm	39.0	2.92	33.2	44.9	30.3	47.8	7.47%	14.94%	22.42%	37.1	41.0
Zr, ppm	9.71	0.725	8.25	11.16	7.53	11.88	7.47%	14.94%	22.42%	9.22	10.19

SI unit equivalents: ppm, parts per million \equiv mg/kg \equiv μ g/g \equiv 0.0001 wt.% \equiv 1000 ppb, parts per billion.

Note 1: intervals may appear asymmetric due to rounding.

Note 2: the number of decimal places quoted does not imply accuracy of the certified value to this level but are given to minimise rounding errors when calculating 2SD and 3SD windows.

Table 2. Indicative Values for OREAS 550.

Constituent	Unit	Value	Constituent	Unit	Value	Constituent	Unit	Value
Borate Fusio	on XRF							
As	ppm	48.3	Мо	ppm	31.5	Sn	ppm	< 10
Bi	ppm	< 10	Na₂O	wt.%	0.048	V ₂ O ₅	ppm	127
Cd	ppm	< 10	Ni	ppm	18.9	Zn	ppm	128
CI	ppm	235	Pb	ppm	< 90	Zr	ppm	83
Cr ₂ O ₃	ppm	69	Sb	ppm	< 10			
Thermograv	imetry							
H ₂ O-	wt.%	0.395						
Peroxide Fu	sion ICP							
As	ppm	8.37	In	ppm	0.22	Та	ppm	0.77
Be	ppm	1.61	Мо	ppm	3.28	Te	ppm	< 1
Cd	ppm	< 10	Ni	ppm	29.8	TI	ppm	< 0.5
Cr	ppm	64	Re	ppm	< 0.1	W	ppm	2.99
Ge	ppm	0.78	Sc	ppm	5.42	Zr	ppm	104
Hf	ppm	3.07	Se	ppm	7.83			
Hg	ppm	0.080	Sn	ppm	2.61			
4-Acid Diges	stion							
В	ppm	11.9	Se	ppm	1.50	Ti	wt.%	0.117
Ge	ppm	0.074	Та	ppm	0.42			
Hg	ppm	0.047	Te	ppm	0.061			
Aqua Regia	Digestio	n						
Au	ppm	0.006	Nb	ppm	0.072	Sm	ppm	2.55
В	ppm	9.96	Pd	ppb	< 10	Та	ppm	< 0.01
Ge	ppm	0.041	Pt	ppb	20.0	Tm	ppm	0.11
Na	wt.%	0.013	Si	wt.%	0.091			

SI unit equivalents: ppm (parts per million) \equiv mg/kg \equiv μ g/g \equiv 0.0001 wt.% \equiv 1000 ppb (parts per billion).

Note: the number of significant figures reported is not a reflection of the level of certainty of stated values. They are instead an artefact of ORE's in-house CRM-specific LIMS.

INTRODUCTION

OREAS reference materials are intended to provide a low-cost method of evaluating and improving the quality of analysis of geological samples. To the geologist they provide a means of implementing quality control in analytical data sets generated in exploration from the grass roots level through to prospect evaluation, and in grade control at mining operations. To the analyst they provide an effective means of calibrating analytical equipment, assessing new techniques and routinely monitoring in-house procedures. OREAS reference materials enable users to successfully achieve process control of these tasks because the observed variance from repeated analysis has its origin almost exclusively in the analytical process rather than the reference material itself. In evaluating laboratory performance with this CRM, the section headed 'Instructions for correct use' should be read carefully.

Table 1 provides performance gate intervals for the certified values. Table 2 shows indicative values, Table 3 provides some indicative physical properties and Table 4 presents the 95% confidence and tolerance limits for all certified values. Tabulated results of all elements together with uncorrected means, medians, standard deviations, relative standard deviations and per cent deviation of lab means from the corrected mean of means (PDM³) are presented in the detailed certification data for this CRM (OREAS 550-DataPack.1.1.210420_162423.xlsx).

Results are also presented in scatter plots for Co and Cu by 4-acid digestion and peroxide fusion ICP in Figures 1 to 4 respectively, together with ±3SD (magenta) and ±5% (yellow) control lines and certified value (green line). Accepted individual results are coloured blue and individual and dataset outliers are identified in red and violet, respectively.

SOURCE MATERIAL

OREAS 550 was prepared from a blend of copper-cobalt oxide ore materials supplied by Metorex (subsidiary of Jinchuan Group) from the Dilala Est Copper Cobalt Project, located within the Kolwezi Klippe area of Lualaba Province (Democratic Republic of the Congo). Kolwezi is ~320 kilometres northwest from Lubumbashi, the provincial capital of Haut Katanga Province. Copper-cobalt ores within the Kolwezi Klippe area are stratiform, sediment-hosted deposits hosted within Neoproterozoic sedimentary rocks of the Central African Copperbelt. Primary mineralisation typically consists of bornite, chalcopyrite and cobalt sulphides hosted within dolomitic argillite, shales, sandstone and dolomite horizons. Significant oxide minerals include malachite, chalcocite, digenite and hetrogenite.

PERFORMANCE GATES

Table 1 above shows intervals calculated for two and three standard deviations. As a guide these intervals may be regarded as warning or rejection for multiple 2SD outliers, or rejection for individual 3SD outliers in QC monitoring, although their precise application should be at the discretion of the QC manager concerned (also see 'Intended Use' section below). Westgard Rules extend the basics of single-rule QC monitoring using multi-rules (for more information visit www.westgard.com/mltirule.htm). A second method utilises a 5% window calculated directly from the certified value.



Standard deviation is also shown in relative percent for one, two and three relative standard deviations (1RSD, 2RSD and 3RSD) to facilitate an appreciation of the magnitude of these numbers and a comparison with the 5% window. Caution should be exercised when concentration levels approach lower limits of detection of the analytical methods employed as performance gates calculated from standard deviations tend to be excessively wide whereas those determined by the 5% method are too narrow. One approach used at commercial laboratories is to set the acceptance criteria at twice the detection level (DL) \pm 10%.

i.e. Certified Value ± 10% ± 2DL (adapted from Govett, 1983).

COMMINUTION AND HOMOGENISATION PROCEDURES

The material constituting OREAS 550 was prepared in the following manner:

- Drying to constant mass at 105° Celsius;
- Multi-stage milling of each source material to achieve a particle size distribution of >98% passing 75 microns;
- Preliminary homogenisation of each source material;
- Representative sampling and check assaying of each source material;
- Blending in appropriate proportions to achieve target grades;
- Packaging in 10g units sealed in laminated foil pouches and 1000g units in plastic iars.

PHYSICAL PROPERTIES

OREAS 550 was tested at ORE Research & Exploration Pty Ltd's onsite facility for various physical properties. Table 3 presents these findings that should be used for informational purposes only.

Table 3. Physical properties of OREAS 550.

Bulk Density (g/L)	Moisture%	Munsell Notation [‡]	Munsell Color‡
850	0.39	5YR 8/1	Pinkish Gray

[‡]The Munsell Rock Color Chart helps geologists and archeologists communicate with colour more effectively by cross-referencing ISCC-NBS colour names with unique Munsell alpha-numeric colour notations for rock colour samples.

ANALYTICAL PROGRAM

Thirty analytical laboratories participated in the program to characterise the elements reported in Table 1. The following methods were employed:

- Lithium borate fusion whole rock analysis package by X-ray fluorescence (up to 19 laboratories depending on the element);
- Thermogravimetry: Loss on Ignition (LOI) at 1000°C (13 laboratories used a thermogravimetric analyser, 7 laboratories included LOI with their fusion package and 4 laboratories used conventional muffle furnace);
- C and S by infrared combustion furnace/CS analyser (27 laboratories);
- Cu and Co by *5% sulphuric acid leach with ICP or AAS finish (up to 19 laboratories);

- Cu and Co by *10% sulphuric acid leach with ICP or AAS finish (up to 12 laboratories);
- Sodium peroxide fusion with full suite elemental package by ICP-OES and/or MS finish (up to 20 laboratories depending on the element);
- 4-acid digestion for full suite elemental package by ICP-OES and MS finish (up to 26 laboratories depending on the element);
- Aqua regia digestion for full elemental suite ICP-OES and ICP-MS (up to 26 laboratories depending on the element).

For the round robin program ten 400g test units were taken at predetermined intervals immediately following homogenisation and are considered representative of the entire prepared batch. The six samples received by each laboratory were obtained by taking two 25g scoop splits from each of three separate 400g test units. This format enabled nested ANOVA treatment of the results to evaluate homogeneity, i.e., to ascertain whether between-unit variance is greater than within-unit variance (see 'Homogeneity Evaluation' section below).

STATISTICAL ANALYSIS

Standard Deviation intervals (see Table 1) provide an indication of a level of performance that might reasonably be expected from a laboratory being monitored by this CRM in a QA/QC program. They take into account errors attributable to measurement uncertainty and CRM variability. For an effective CRM the contribution of the latter should be negligible in comparison to measurement errors. The Standard Deviation values include all sources of measurement uncertainty: between-laboratory variance, within-run variance (precision errors) and CRM variability.

In the application of SD's in monitoring performance it is important to note that not all laboratories function at the same level of proficiency and that different methods in use at a particular laboratory have differing levels of precision. Each laboratory has its own inherent SD (for a specific concentration level and analyte-method pair) based on the analytical process and this SD is not directly related to the round robin program (see Intended Use section for more detail).

The SD for each analyte's certified value is calculated from the same filtered data set used to determine the certified value, i.e., after removal of all individual, lab dataset (batch) and 3SD outliers (single iteration). These outliers can only be removed after the absolute homogeneity of the CRM has been independently established, i.e., the outliers must be confidently deemed to be analytical rather than arising from inhomogeneity of the CRM. The standard deviation is then calculated for each analyte from the pooled accepted analyses generated from the certification program.

Certified Values, Confidence Limits and Tolerance Limits (Table 4) have been determined for each analyte following removal of individual, laboratory dataset (batch) and 3SD outliers (single iteration).

95% Confidence Limits are inversely proportional to the number of participating laboratories and inter-laboratory agreement. It is a measure of the reliability of the certified

^{*}See 'Appendix' for specified methodology.

value. A 95% confidence interval indicates a 95% probability that the true value of the analyte under consideration lies between the upper and lower limits. **95% Confidence** *Limits should not be used as control limits for laboratory performance.*

Indicative (uncertified) values (Table 2) are provided where i) a laboratory reported analytes beyond those requested; ii) the number of laboratories reporting a particular analyte is insufficient (< 5) to support certification; iii) inter-laboratory consensus is poor; or iv) a significant proportion of results are outlying.

For individual outliers within a laboratory batch the z-score test is used in combination with a second method that determines the per cent deviation of the individual value from the batch median. Outliers in general are selected on the basis of z-scores > 2.5 and with per cent deviations (i) > 3 and (ii) more than three times the average absolute per cent deviation for the batch. In certain instances, statistician's prerogative has been employed in discriminating outliers. Each laboratory data set mean is tested for outlying status based on z-score discrimination and rejected if > 2.5. After individual and laboratory data set (batch) outliers have been eliminated a non-iterative 3 standard deviation filter is applied, with those values lying outside this window also relegated to outlying status.

Table 4. 95% Confidence & Tolerance Limits for OREAS 550.

Constituent	Certified	95% Confid	ence Limits	95% Toler	ance Limits
Constituent	Value	Low	High	Low	High
Borate Fusion XRF					
Al ₂ O ₃ , Aluminium(III) oxide (wt.%)	5.39	5.36	5.42	5.35	5.43
BaO, Barium oxide (ppm)	386	364	407	343	429
CaO, Calcium oxide (wt.%)	9.70	9.64	9.75	9.63	9.76
Co, Cobalt (wt.%)	0.145	0.143	0.147	0.143	0.147
Cu, Copper (wt.%)	0.278	0.273	0.284	0.271	0.285
Fe ₂ O ₃ , Iron(III) oxide (wt.%)	1.89	1.88	1.90	1.88	1.91
K ₂ O, Potassium oxide (wt.%)	0.709	0.703	0.715	0.700	0.718
MgO, Magnesium oxide (wt.%)	10.72	10.67	10.77	10.64	10.81
MnO, Manganese oxide (wt.%)	0.160	0.159	0.161	0.158	0.162
P ₂ O ₅ , Phosphorus(V) oxide (wt.%)	0.081	0.079	0.084	0.080	0.083
SiO ₂ , Silicon dioxide (wt.%)	52.79	52.66	52.93	52.58	53.01
SO ₃ , Sulphur trioxide (wt.%)	0.565	0.554	0.576	0.552	0.578
SrO, Strontium oxide (ppm)	< 100	IND	IND	IND	IND
TiO ₂ , Titanium dioxide (wt.%)	0.314	0.310	0.317	0.305	0.322
Thermogravimetry					
LOI ¹⁰⁰⁰ , Loss On Ignition @1000°C (wt.%)	16.91	16.85	16.97	16.84	16.97
Infrared Combustion					
C, Carbon (wt.%)	4.11	4.07	4.15	4.08	4.14
S, Sulphur (wt.%)	0.220	0.215	0.225	0.212	0.228
Sulphuric Acid 5% Leach					
Co, Cobalt (wt.%)	0.101	0.088	0.115	0.099	0.104
Cu, Copper (wt.%)	0.197	0.191	0.203	0.193	0.201

SI unit equivalents: ppm, parts per million \equiv mg/kg \equiv μ g/g \equiv 0.0001 wt.% \equiv 1000 ppb, parts per billion. Note: intervals may appear asymmetric due to rounding; IND = indeterminate (due to reading resolution error).

Table 4 continued. Certified 95% Confidence Limits 95% Tolerance Limits										
Constituent	Value	Low	High	Low	High					
Sulphuric Acid 10% Leach	value	LOW	riigii	LOW	Tilgii					
Co, Cobalt (wt.%)	0.097	0.087	0.107	0.094	0.100					
Cu, Copper (wt.%)	0.201	0.186	0.217	0.196	0.207					
Peroxide Fusion ICP	0.201	0.100	0.211	0.100	0.207					
Ag, Silver (ppm)	< 1	IND	IND	IND	IND					
Al, Aluminium (wt.%)	2.87	2.85	2.89	2.82	2.92					
B, Boron (ppm)	229	219	238	219	239					
Ba, Barium (ppm)	349	343	354	336	362					
Bi, Bismuth (ppm)	1.45	1.24	1.66	IND	IND					
Ca, Calcium (wt.%)	7.03	6.97	7.08	6.89	7.16					
Ce, Cerium (ppm)	39.3	37.6	41.1	37.5	41.2					
Co, Cobalt (wt.%)	0.150	0.148	0.152	0.147	0.153					
Cs, Caesium (ppm)	0.60	0.57	0.64	IND	IND					
Cu, Copper (wt.%)	0.279	0.276	0.283	0.272	0.286					
Dy, Dysprosium (ppm)	3.57	3.39	3.76	3.31	3.84					
Er, Erbium (ppm)	1.93	1.80	2.06	1.81	2.05					
Eu, Europium (ppm)	0.92	0.87	0.96	0.82	1.02					
Fe, Iron (wt.%)	1.33	1.31	1.34	1.29	1.36					
Ga, Gallium (ppm)	7.80	7.36	8.25	7.40	8.20					
Gd, Gadolinium (ppm)	4.14	3.94	4.34	3.91	4.37					
Ho, Holmium (ppm)	0.66	0.59	0.73	0.61	0.70					
K, Potassium (wt.%)	0.628	0.610	0.645	0.602	0.653					
La, Lanthanum (ppm)	22.8	21.1	24.5	21.8	23.9					
Li, Lithium (ppm)	105	100	110	101	109					
Lu, Lutetium (ppm)	0.26	0.23	0.28	0.22	0.29					
Mg, Magnesium (wt.%)	6.44	6.41	6.47	6.34	6.53					
Mn, Manganese (wt.%)	0.123	0.121	0.125	0.119	0.126					
Nb, Niobium (ppm)	9.37	8.36	10.38	IND	IND					
Nd, Neodymium (ppm)	19.4	17.8	20.9	18.2	20.5					
P, Phosphorus (wt.%)	0.038	0.033	0.043	IND	IND					
Pb, Lead (ppm)	16.8	13.0	20.7	14.3	19.3					
Pr, Praseodymium (ppm)	5.13	4.83	5.42	4.92	5.33					
Rb, Rubidium (ppm)	22.9	22.1	23.8	21.7	24.1					
S, Sulphur (wt.%)	0.221	0.210	0.232	0.205	0.237					
Sb, Antimony (ppm)	< 1	IND	IND	IND	IND					
Si, Silicon (wt.%)	25.30	24.94	25.66	24.85	25.75					
Sm, Samarium (ppm)	4.37	4.17	4.56	4.01	4.72					
Sr, Strontium (ppm)	59	57	61	57	61					
Tb, Terbium (ppm)	0.61	0.57	0.64	0.57	0.64					

SI unit equivalents: ppm, parts per million \equiv mg/kg \equiv μ g/g \equiv 0.0001 wt.% \equiv 1000 ppb, parts per billion.

Note: intervals may appear asymmetric due to rounding; IND = indeterminate (due to reading resolution error).



	Table 4 contin	1			
Constituent	Certified		lence Limits		ance Limits
	Value	Low	High	Low	High
Peroxide Fusion ICP continued		T			T
Th, Thorium (ppm)	5.37	5.11	5.62	5.10	5.63
Ti, Titanium (wt.%)	0.188	0.184	0.191	0.180	0.195
Tm, Thulium (ppm)	0.27	0.25	0.29	0.25	0.29
U, Uranium (ppm)	23.9	22.3	25.5	22.4	25.4
V, Vanadium (ppm)	56	51	61	52	61
Y, Yttrium (ppm)	17.8	17.0	18.6	17.1	18.5
Yb, Ytterbium (ppm)	1.70	1.60	1.81	IND	IND
Zn, Zinc (ppm)	51	46	57	45	58
*4-Acid Digestion					
Ag, Silver (ppm)	0.127	0.119	0.136	IND	IND
Al, Aluminium (wt.%)	2.84	2.81	2.86	2.78	2.89
As, Arsenic (ppm)	5.82	5.51	6.12	5.40	6.23
Ba, Barium (ppm)	352	348	356	343	360
Be, Beryllium (ppm)	1.31	1.27	1.35	1.27	1.35
Bi, Bismuth (ppm)	1.40	1.37	1.44	1.37	1.44
Ca, Calcium (wt.%)	6.91	6.82	7.00	6.80	7.02
Cd, Cadmium (ppm)	0.13	0.13	0.14	IND	IND
Ce, Cerium (ppm)	36.6	35.1	38.1	35.4	37.7
Co, Cobalt (wt.%)	0.148	0.144	0.151	0.145	0.151
Cr, Chromium (ppm)	23.1	20.2	25.9	21.5	24.6
Cs, Caesium (ppm)	0.58	0.56	0.60	0.55	0.61
Cu, Copper (wt.%)	0.277	0.273	0.281	0.271	0.283
Dy, Dysprosium (ppm)	3.30	3.24	3.36	3.18	3.43
Er, Erbium (ppm)	1.76	1.68	1.84	1.67	1.85
Eu, Europium (ppm)	0.84	0.82	0.85	0.80	0.87
Fe, Iron (wt.%)	1.28	1.25	1.30	1.25	1.31
Ga, Gallium (ppm)	7.65	7.53	7.77	7.42	7.87
Gd, Gadolinium (ppm)	4.02	3.86	4.18	3.81	4.24
Hf, Hafnium (ppm)	2.19	2.10	2.28	2.08	2.30
Ho, Holmium (ppm)	0.62	0.59	0.66	0.60	0.64
In, Indium (ppm)	0.18	0.17	0.18	0.17	0.18
K, Potassium (wt.%)	0.614	0.599	0.628	0.600	0.628
La, Lanthanum (ppm)	19.2	18.4	20.1	18.7	19.8
Li, Lithium (ppm)	105	102	108	102	108
Lu, Lutetium (ppm)	0.23	0.22	0.25	0.22	0.25
Mg, Magnesium (wt.%)	6.43	6.33	6.52	6.29	6.56
Mn, Manganese (wt.%)	0.119	0.117	0.122	0.118	0.121
L	_ t	-1			I

SI unit equivalents: ppm, parts per million \equiv mg/kg \equiv μ g/g \equiv 0.0001 wt.% \equiv 1000 ppb, parts per billion.

Note: intervals may appear asymmetric due to rounding; IND = indeterminate (due to reading resolution error).



^{*}Four acid digestion quantitatively dissolves nearly all minerals in the majority of geological samples however, some refractory minerals may only be partially digested.

Table 4 continued.										
Constituent	Certified		dence Limits		ance Limits					
	Value	Low	High	Low	High					
*4-Acid Digestion continued	T				I					
Mo, Molybdenum (ppm)	2.25	2.21	2.30	2.16	2.35					
Na, Sodium (wt.%)	0.040	0.036	0.044	0.038	0.042					
Nb, Niobium (ppm)	5.98	5.37	6.58	5.72	6.23					
Nd, Neodymium (ppm)	18.5	17.9	19.0	18.0	18.9					
Ni, Nickel (ppm)	16.7	16.1	17.2	16.1	17.3					
P, Phosphorus (wt.%)	0.036	0.035	0.037	0.035	0.038					
Pb, Lead (ppm)	13.7	13.2	14.3	13.0	14.5					
Pr, Praseodymium (ppm)	4.61	4.40	4.82	4.47	4.75					
Rb, Rubidium (ppm)	22.7	22.1	23.3	22.1	23.3					
Re, Rhenium (ppm)	0.008	0.007	0.008	IND	IND					
S, Sulphur (wt.%)	0.233	0.228	0.239	0.224	0.243					
Sb, Antimony (ppm)	0.32	0.30	0.33	0.29	0.35					
Sc, Scandium (ppm)	6.27	6.11	6.43	6.10	6.43					
Sm, Samarium (ppm)	4.09	3.88	4.29	3.85	4.33					
Sn, Tin (ppm)	1.18	1.13	1.24	IND	IND					
Sr, Strontium (ppm)	56	55	58	55	58					
Tb, Terbium (ppm)	0.62	0.60	0.65	0.60	0.64					
Th, Thorium (ppm)	5.15	4.96	5.34	4.99	5.31					
TI, Thallium (ppm)	0.094	0.091	0.097	IND	IND					
Tm, Thulium (ppm)	0.25	0.23	0.27	0.24	0.27					
U, Uranium (ppm)	24.7	23.7	25.7	23.1	26.3					
V, Vanadium (ppm)	53	51	55	51	55					
W, Tungsten (ppm)	0.78	0.73	0.84	IND	IND					
Y, Yttrium (ppm)	16.5	16.0	17.0	16.1	16.9					
Yb, Ytterbium (ppm)	1.62	1.57	1.67	1.52	1.72					
Zn, Zinc (ppm)	41.9	40.8	43.0	39.8	44.1					
Zr, Zirconium (ppm)	76	74	78	74	78					
Aqua Regia Digestion		_								
Ag, Silver (ppm)	0.086	0.082	0.091	0.070	0.102					
Al, Aluminium (wt.%)	1.55	1.50	1.59	1.50	1.59					
As, Arsenic (ppm)	4.93	4.75	5.11	4.68	5.18					
Ba, Barium (ppm)	297	293	300	288	305					
Be, Beryllium (ppm)	0.92	0.90	0.95	0.89	0.96					
Bi, Bismuth (ppm)	1.07	1.05	1.09	1.03	1.11					
Ca, Calcium (wt.%)	6.70	6.62	6.79	6.60	6.80					
Cd, Cadmium (ppm)	0.13	0.13	0.13	IND	IND					
Ce, Cerium (ppm)	15.8	15.2	16.4	15.3	16.3					
/	1	1	1		L					

SI unit equivalents: ppm, parts per million \equiv mg/kg \equiv μ g/g \equiv 0.0001 wt.% \equiv 1000 ppb, parts per billion.

Note: intervals may appear asymmetric due to rounding; IND = indeterminate (due to reading resolution error).



^{*}Four acid digestion quantitatively dissolves nearly all minerals in the majority of geological samples however, some refractory minerals may only be partially digested.

Constituent Certified Value β5% Coloration Coloration (Date of Value) Low (Date of Value) Low (Date of Value) Home (Date of Value) Ho	Table 4 continued.							
Aqua Regia Digestion continued Low High Low Lind Lind Col. Low Low Hind Low Low Hind Low Low Low Low Low <td>Constituent</td> <td></td> <td colspan="2">95% Confidence Limits</td> <td colspan="2">95% Tolerance Limits</td>	Constituent		95% Confidence Limits		95% Tolerance Limits			
Co, Cobalt (wt.%) 0.139 0.136 0.143 0.136 0.142 Cr, Chromium (ppm) 21.5 20.9 22.0 20.4 22.6 Cs, Caesium (ppm) 0.13 0.11 0.15 IND IND Cu, Copper (wt.%) 0.274 0.271 0.277 0.269 0.280 Dy, Dysprosium (ppm) 1.92 1.84 2.01 IND IND Er, Erbium (ppm) 0.85 0.81 0.90 IND IND Eu, Europium (ppm) 0.56 0.53 0.58 IND IND Eu, Europium (ppm) 1.07 1.05 1.10 1.05 1.10 Ga, Gallium (ppm) 4.45 4.30 4.61 4.31 4.60 Gd, Gadolinium (ppm) 2.70 2.34 3.06 2.50 2.91 Hf, Hafnium (ppm) 0.31 0.28 0.34 0.29 0.32 Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND In, Indium (ppm) 0.16	Constituent	Value	Low	High	Low	High		
Cr, Chromium (ppm) 21.5 20.9 22.0 20.4 22.6 Cs, Caesium (ppm) 0.13 0.11 0.15 IND IND Cu, Copper (wt.%) 0.274 0.271 0.277 0.269 0.280 Dy, Dysprosium (ppm) 1.92 1.84 2.01 IND IND Eu, Europium (ppm) 0.85 0.81 0.90 IND IND Eu, Europium (ppm) 0.56 0.53 0.58 IND IND Fe, Iron (wt.%) 1.07 1.05 1.10 1.05 1.10 Ga, Gallium (ppm) 4.45 4.30 4.61 4.31 4.60 Gd, Gadolinium (ppm) 2.70 2.34 3.06 2.50 2.91 Hf, Hafnium (ppm) 0.031 0.28 0.34 0.29 0.32 Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND Ho, Holmium (ppm) 0.046 0.15 0.16 0.15 0.16 K, Potassium (wt.%) 0.126						T		
Cs, Caesium (ppm) 0.13 0.11 0.15 IND IND Cu, Copper (wt.%) 0.274 0.271 0.277 0.269 0.280 Dy, Dysprosium (ppm) 1.92 1.84 2.01 IND IND Er, Erbium (ppm) 0.85 0.81 0.90 IND IND Er, Erbium (ppm) 0.56 0.53 0.58 IND IND Eu, Europium (ppm) 0.56 0.53 0.58 IND IND Fe, Iron (wt.%) 1.07 1.05 1.10 1.05 1.10 Ga, Gallium (ppm) 4.45 4.30 4.61 4.31 4.60 Gd, Gadolinium (ppm) 0.27 2.34 3.06 2.50 2.91 Hf, Hafnium (ppm) 0.31 0.28 0.34 0.29 0.32 Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND In, Indium (ppm) 0.34 0.33 0.35 IND IND In, Hg, Herbinim (ppm) 0.16 <	Co, Cobalt (wt.%)	0.139	0.136	0.143	0.136	0.142		
Cu, Copper (wf.%) 0.274 0.271 0.277 0.269 0.280 Dy, Dysprosium (ppm) 1.92 1.84 2.01 IND IND Er, Erbium (ppm) 0.85 0.81 0.90 IND IND Eu, Europium (ppm) 0.56 0.53 0.58 IND IND Eu, Europium (ppm) 0.56 0.53 0.58 IND IND Eu, Europium (ppm) 0.56 0.53 0.58 IND IND Ga, Gallium (ppm) 4.45 4.30 4.61 4.31 4.60 Gd, Gadolinium (ppm) 2.70 2.34 3.06 2.50 2.91 Hf, Hafnium (ppm) 0.31 0.28 0.34 0.29 0.32 Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND In, Indium (ppm) 0.16 0.15 0.16 0.15 0.15 In, Indium (ppm) 0.16 0.15 0.16 0.15 0.18 In, Lithium (ppm) 0.16	Cr, Chromium (ppm)	21.5	20.9	22.0	20.4	22.6		
Dy, Dysprosium (ppm) 1.92 1.84 2.01 IND IND Er, Erbium (ppm) 0.85 0.81 0.90 IND IND Eu, Europium (ppm) 0.56 0.53 0.58 IND IND Fe, Iron (wt.%) 1.07 1.05 1.10 1.05 1.10 Ga, Gallium (ppm) 4.45 4.30 4.61 4.31 4.60 Gd, Gadolinium (ppm) 2.70 2.34 3.06 2.50 2.91 Hf, Hafnium (ppm) 0.31 0.28 0.34 0.29 0.32 Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND Ho, Holmium (ppm) 0.34 0.33 0.35 IND IND In, Indium (ppm) 0.16 0.15 0.16 0.15 0.16 In, Indium (ppm) 0.16 0.15 0.16 0.15 0.16 Li, Lithium (ppm) 5.31 4.94 5.68 5.14 5.49 Li, Lithium (ppm) 0.099 0.	Cs, Caesium (ppm)	0.13	0.11	0.15	IND	IND		
Er, Erbium (ppm) 0.85 0.81 0.90 IND IND Eu, Europium (ppm) 0.56 0.53 0.58 IND IND Fe, Iron (wt.%) 1.07 1.05 1.10 1.05 1.10 Ga, Gallium (ppm) 4.45 4.30 4.61 4.31 4.60 Gd, Gadolinium (ppm) 2.70 2.34 3.06 2.50 2.91 Hf, Hafnium (ppm) 0.31 0.28 0.34 0.29 0.32 Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND Ho, Holmium (ppm) 0.34 0.33 0.35 IND IND In, Indium (ppm) 0.16 0.15 0.16 0.15 0.16 K, Potassium (wt.%) 0.126 0.119 0.133 0.121 0.132 L, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.018 5.99	Cu, Copper (wt.%)	0.274	0.271	0.277	0.269	0.280		
Eu, Europium (ppm) 0.56 0.53 0.58 IND IND Fe, Iron (wt.%) 1.07 1.05 1.10 1.05 1.10 Ga, Gallium (ppm) 4.45 4.30 4.61 4.31 4.60 Gd, Gadolinium (ppm) 2.70 2.34 3.06 2.50 2.91 Hf, Hafnium (ppm) 0.31 0.28 0.34 0.29 0.32 Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND Ho, Holmium (ppm) 0.34 0.33 0.35 IND IND Ho, Holmium (ppm) 0.16 0.15 0.16 0.15 0.16 K, Potassium (wt.%) 0.126 0.119 0.133 0.121 0.132 La, Lanthanum (ppm) 5.31 4.94 5.68 5.14 5.49 Li, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08	Dy, Dysprosium (ppm)	1.92	1.84	2.01	IND	IND		
Fe, Iron (wt.%) 1.07 1.05 1.10 1.05 1.10 Ga, Gallium (ppm) 4.45 4.30 4.61 4.31 4.60 Gd, Gadolinium (ppm) 2.70 2.34 3.06 2.50 2.91 Hf, Hafnium (ppm) 0.31 0.28 0.34 0.29 0.32 Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND Ho, Holmium (ppm) 0.34 0.33 0.35 IND IND In, Indium (ppm) 0.16 0.15 0.16 0.15 0.16 K, Potassium (wt.%) 0.126 0.119 0.133 0.121 0.132 La, Lanthanum (ppm) 5.31 4.94 5.68 5.14 5.49 Li, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mn, Manganese (wt.%) 0.118	Er, Erbium (ppm)	0.85	0.81	0.90	IND	IND		
Ga, Gallium (ppm) 4.45 4.30 4.61 4.31 4.60 Gd, Gadolinium (ppm) 2.70 2.34 3.06 2.50 2.91 Hf, Hafnium (ppm) 0.31 0.28 0.34 0.29 0.32 Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND Ho, Holmium (ppm) 0.14 0.33 0.35 IND IND Ho, Holmium (ppm) 0.16 0.15 0.16 0.15 0.16 K, Potassium (wt.%) 0.126 0.119 0.133 0.121 0.132 La, Lanthanum (ppm) 86 82 90 84 88 Li, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mn, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.122 Mo, Nolybdenum (ppm) 2.12	Eu, Europium (ppm)	0.56	0.53	0.58	IND	IND		
Gd, Gadolinium (ppm) 2.70 2.34 3.06 2.50 2.91 Hf, Hafnium (ppm) 0.31 0.28 0.34 0.29 0.32 Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND Ho, Holmium (ppm) 0.34 0.33 0.35 IND IND In, Indium (ppm) 0.16 0.15 0.16 0.15 0.16 K, Potassium (wt.%) 0.126 0.119 0.133 0.121 0.132 La, Lanthanum (ppm) 86 82 90 84 88 Li, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mh, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.122 Md, Neodymium (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 15.4	Fe, Iron (wt.%)	1.07	1.05	1.10	1.05	1.10		
Hf, Hafnium (ppm) 0.31 0.28 0.34 0.29 0.32 Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND Ho, Holmium (ppm) 0.34 0.33 0.35 IND IND In, Indium (ppm) 0.16 0.15 0.16 0.15 0.16 K, Potassium (wt.%) 0.126 0.119 0.133 0.121 0.132 La, Lanthanum (ppm) 5.31 4.94 5.68 5.14 5.49 Li, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mn, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.121 Mo, Molybdenum (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 </td <td>Ga, Gallium (ppm)</td> <td>4.45</td> <td>4.30</td> <td>4.61</td> <td>4.31</td> <td>4.60</td>	Ga, Gallium (ppm)	4.45	4.30	4.61	4.31	4.60		
Hg, Mercury (ppm) 0.022 0.018 0.027 IND IND Ho, Holmium (ppm) 0.34 0.33 0.35 IND IND In, Indium (ppm) 0.16 0.15 0.16 0.15 0.16 K, Potassium (wt.%) 0.126 0.119 0.133 0.121 0.132 La, Lanthanum (ppm) 5.31 4.94 5.68 5.14 5.49 Li, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mn, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.121 Mo, Molybdenum (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.0	Gd, Gadolinium (ppm)	2.70	2.34	3.06	2.50	2.91		
Ho, Holmium (ppm) 0.34 0.33 0.35 IND IND In, Indium (ppm) 0.16 0.15 0.16 0.15 0.16 K, Potassium (wt.%) 0.126 0.119 0.133 0.121 0.132 La, Lanthanum (ppm) 5.31 4.94 5.68 5.14 5.49 Li, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mn, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.121 Mo, Molybdenum (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11	Hf, Hafnium (ppm)	0.31	0.28	0.34	0.29	0.32		
In, Indium (ppm) 0.16 0.15 0.16 0.15 0.16 K, Potassium (wt.%) 0.126 0.119 0.133 0.121 0.132 La, Lanthanum (ppm) 5.31 4.94 5.68 5.14 5.49 Li, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mn, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.121 Mo, Molybdenum (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm)	Hg, Mercury (ppm)	0.022	0.018	0.027	IND	IND		
K, Potassium (wt.%) 0.126 0.119 0.133 0.121 0.132 La, Lanthanum (ppm) 5.31 4.94 5.68 5.14 5.49 Li, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mn, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.121 Mo, Molybdenum (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Reneium (p	Ho, Holmium (ppm)	0.34	0.33	0.35	IND	IND		
La, Lanthanum (ppm) 5.31 4.94 5.68 5.14 5.49 Li, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mn, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.121 Mo, Molybdenum (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) <td>In, Indium (ppm)</td> <td>0.16</td> <td>0.15</td> <td>0.16</td> <td>0.15</td> <td>0.16</td>	In, Indium (ppm)	0.16	0.15	0.16	0.15	0.16		
Li, Lithium (ppm) 86 82 90 84 88 Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mn, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.121 Mo, Molybdenum (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.2	K, Potassium (wt.%)	0.126	0.119	0.133	0.121	0.132		
Lu, Lutetium (ppm) 0.099 0.090 0.108 IND IND Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mn, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.121 Mo, Molybdenum (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimory (ppm)	La, Lanthanum (ppm)	5.31	4.94	5.68	5.14	5.49		
Mg, Magnesium (wt.%) 6.08 5.99 6.18 5.97 6.20 Mn, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.121 Mo, Molybdenum (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm)	Li, Lithium (ppm)	86	82	90	84	88		
Mn, Manganese (wt.%) 0.118 0.115 0.122 0.115 0.121 Mo, Molybdenum (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm)	Lu, Lutetium (ppm)	0.099	0.090	0.108	IND	IND		
Mo, Molybdenum (ppm) 2.12 2.07 2.17 2.04 2.20 Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm)	Mg, Magnesium (wt.%)	6.08	5.99	6.18	5.97	6.20		
Nd, Neodymium (ppm) 8.32 7.70 8.93 7.50 9.13 Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm) 1.11 0.98 1.25 IND IND Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm)	Mn, Manganese (wt.%)	0.118	0.115	0.122	0.115	0.121		
Ni, Nickel (ppm) 15.4 14.8 16.0 14.9 15.8 P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm) 1.11 0.98 1.25 IND IND Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm)	Mo, Molybdenum (ppm)	2.12	2.07	2.17	2.04	2.20		
P, Phosphorus (wt.%) 0.033 0.032 0.034 0.032 0.034 Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm) 1.11 0.98 1.25 IND IND Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23<	Nd, Neodymium (ppm)	8.32	7.70	8.93	7.50	9.13		
Pb, Lead (ppm) 11.9 11.4 12.5 11.4 12.5 Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm) 1.11 0.98 1.25 IND IND Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33 <td>Ni, Nickel (ppm)</td> <td>15.4</td> <td>14.8</td> <td>16.0</td> <td>14.9</td> <td>15.8</td>	Ni, Nickel (ppm)	15.4	14.8	16.0	14.9	15.8		
Pr, Praseodymium (ppm) 1.60 1.25 1.96 IND IND Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm) 1.11 0.98 1.25 IND IND Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	P, Phosphorus (wt.%)	0.033	0.032	0.034	0.032	0.034		
Rb, Rubidium (ppm) 4.16 3.89 4.42 3.98 4.34 Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm) 1.11 0.98 1.25 IND IND Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.035 0.031 0.040 0.36 0.41 Te, Tellurium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	Pb, Lead (ppm)	11.9	11.4	12.5	11.4	12.5		
Re, Rhenium (ppm) 0.007 0.007 0.008 IND IND S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm) 1.11 0.98 1.25 IND IND Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.035 0.031 0.040 0.36 0.41 Te, Tellurium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	Pr, Praseodymium (ppm)	1.60	1.25	1.96	IND	IND		
S, Sulphur (wt.%) 0.224 0.215 0.232 0.216 0.231 Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm) 1.11 0.98 1.25 IND IND Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.39 0.38 0.40 0.36 0.41 Te, Tellurium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	Rb, Rubidium (ppm)	4.16	3.89	4.42	3.98	4.34		
Sb, Antimony (ppm) 0.19 0.17 0.20 0.17 0.21 Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm) 1.11 0.98 1.25 IND IND Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.39 0.38 0.40 0.36 0.41 Te, Tellurium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	Re, Rhenium (ppm)	0.007	0.007	0.008	IND	IND		
Sc, Scandium (ppm) 4.44 4.25 4.63 4.28 4.60 Se, Selenium (ppm) 1.11 0.98 1.25 IND IND Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.39 0.38 0.40 0.36 0.41 Te, Tellurium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	S, Sulphur (wt.%)	0.224	0.215	0.232	0.216	0.231		
Se, Selenium (ppm) 1.11 0.98 1.25 IND IND Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.39 0.38 0.40 0.36 0.41 Te, Tellurium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	Sb, Antimony (ppm)	0.19	0.17	0.20	0.17	0.21		
Sn, Tin (ppm) 0.35 0.31 0.39 0.31 0.40 Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.39 0.38 0.40 0.36 0.41 Te, Tellurium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	Sc, Scandium (ppm)	4.44	4.25	4.63	4.28	4.60		
Sr, Strontium (ppm) 44.3 42.8 45.7 43.4 45.1 Tb, Terbium (ppm) 0.39 0.38 0.40 0.36 0.41 Te, Tellurium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	Se, Selenium (ppm)	1.11	0.98	1.25	IND	IND		
Tb, Terbium (ppm) 0.39 0.38 0.40 0.36 0.41 Te, Tellurium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	Sn, Tin (ppm)	0.35	0.31	0.39	0.31	0.40		
Te, Tellurium (ppm) 0.035 0.031 0.040 IND IND Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	Sr, Strontium (ppm)	44.3	42.8	45.7	43.4	45.1		
Th, Thorium (ppm) 2.23 2.15 2.30 2.12 2.33	Tb, Terbium (ppm)	0.39	0.38	0.40	0.36	0.41		
	Te, Tellurium (ppm)	0.035	0.031	0.040	IND	IND		
Ti, Titanium (wt.%) < 0.01 IND IND IND	Th, Thorium (ppm)	2.23	2.15	2.30	2.12	2.33		
	Ti, Titanium (wt.%)	< 0.01	IND	IND	IND	IND		

SI unit equivalents: ppm, parts per million \equiv mg/kg \equiv µg/g \equiv 0.0001 wt.% \equiv 1000 ppb, parts per billion. Note: intervals may appear asymmetric due to rounding; IND = indeterminate (due to reading resolution error).



Constituent	Certified	95% Confidence Limits		95% Tolerance Limits				
	Value	Low	High	Low	High			
Aqua Regia Digestion continued								
TI, Thallium (ppm)	0.050	0.048	0.051	IND	IND			
U, Uranium (ppm)	22.3	21.8	22.8	20.8	23.7			
V, Vanadium (ppm)	33.1	32.3	33.8	32.1	34.1			
W, Tungsten (ppm)	0.28	0.26	0.30	0.23	0.33			
Y, Yttrium (ppm)	8.81	8.58	9.04	8.62	9.00			
Yb, Ytterbium (ppm)	0.71	0.68	0.74	IND	IND			
Zn, Zinc (ppm)	39.0	37.8	40.3	37.6	40.5			
Zr, Zirconium (ppm)	9.71	9.29	10.12	9.37	10.04			

SI unit equivalents: ppm, parts per million \equiv mg/kg \equiv μ g/g \equiv 0.0001 wt.% \equiv 1000 ppb, parts per billion. Note: intervals may appear asymmetric due to rounding; IND = indeterminate (due to reading resolution error).

Homogeneity Evaluation

The tolerance limits (ISO 16269:2014) shown in Table 4 were determined using an analysis of precision errors method and are considered a conservative estimate of true homogeneity. The meaning of tolerance limits may be illustrated for copper by fusion with XRF, where 99% of the time $(1-\alpha=0.99)$ at least 95% of subsamples $(\rho=0.95)$ will have concentrations lying between 0.271 and 0.285 wt.%. Put more precisely, this means that if the same number of subsamples were taken and analysed in the same manner repeatedly, 99% of the tolerance intervals so constructed would cover at least 95% of the total population, and 1% of the tolerance intervals would cover less than 95% of the total population (ISO Guide 35). *Please note that tolerance limits pertain to the homogeneity of the CRM only and should not be used as control limits for laboratory performance.*

Based on the statistical analysis of the results of the inter-laboratory certification program it can be concluded that OREAS 550 is sufficiently homogenous and fit-for-purpose as a certified reference material (see 'Intended Use' below).

PARTICIPATING LABORATORIES

- 1. Actlabs, Ancaster, Ontario, Canada
- 2. AGAT Laboratories, Mississauga, Ontario, Canada
- 3. AH Knight, Kitwe, Copperbelt, Zambia
- 4. AH Knight, Lumwana, Mwinilunga, Zambia
- 5. ALS, Brisbane, QLD, Australia
- 6. ALS, Johannesburg, South Africa
- 7. ALS, Lima, Peru
- 8. ALS, Loughrea, Galway, Ireland
- 9. ALS, Vancouver, BC, Canada
- 10. American Assay Laboratories, Sparks, Nevada, USA

- 11. Bureau Veritas Commodities Canada Ltd, Vancouver, BC, Canada
- 12. Bureau Veritas Geoanalytical, Adelaide, SA, Australia
- 13. Bureau Veritas Geoanalytical, Perth, WA, Australia
- 14. Carsurin, Kendari, Sulawesi, Indonesia
- 15. Inspectorate (BV), Lima, Peru
- 16. Intertek Genalysis, Perth, WA, Australia
- 17. Intertek Testing Services Philippines, Cupang, Muntinlupa, Philippines
- 18. MSALABS, Vancouver, BC, Canada
- 19. Nagrom, Perth, WA, Australia
- 20. Ontario Geological Survey, Sudbury, Ontario, Canada
- 21. PT Geoservices Ltd, Cikarang, Jakarta Raya, Indonesia
- 22. PT Intertek Utama Services, Jakarta Timur, DKI Jakarta, Indonesia
- 23. Robinson International Afrique SAS, Lubumbashi, Katanga, Democratic Republic of the Congo
- 24. Ruashi Mine Site Laboratory, Lubumbashi, Katanga, Democratic Republic of the Congo
- 25. Saskatchewan Research Council, Saskatoon, Saskatchewan, Canada
- 26. SGS, Randfontein, Gauteng, South Africa
- 27. SGS Canada Inc., Vancouver, BC, Canada
- 28. SGS del Peru, Lima, Peru
- 29. SGS Geosol Laboratorios Ltda, Vespasiano, Minas Gerais, Brazil
- 30. UIS Analytical Services, Centurion, South Africa

Please note: To preserve anonymity, the above numbered alphabetical list of participating laboratories does not correspond with the Lab ID numbering on the scatter plots below.

PREPARER AND SUPPLIER

Certified reference material OREAS 550 was prepared and certified by:




ORE Research & Exploration Pty Ltd


Tel: +613-9729 0333
37A Hosie Street


Fax: +613-9729 8338
Bayswater North VIC 3153


Web: www.oreas.com
AUSTRALIA

Email: info@ore.com.au

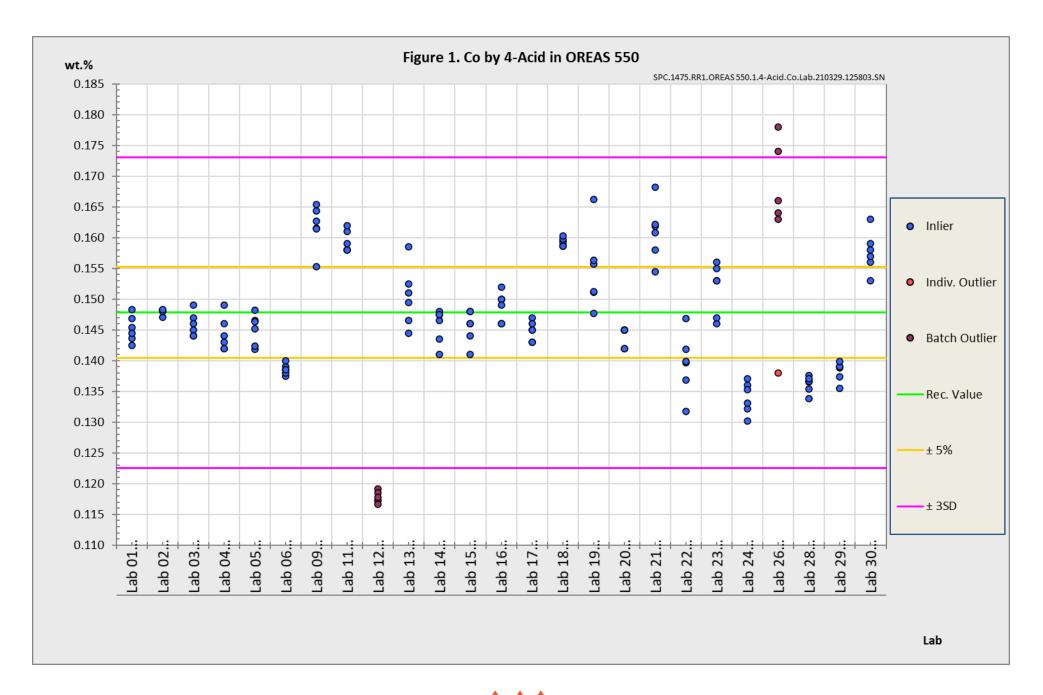

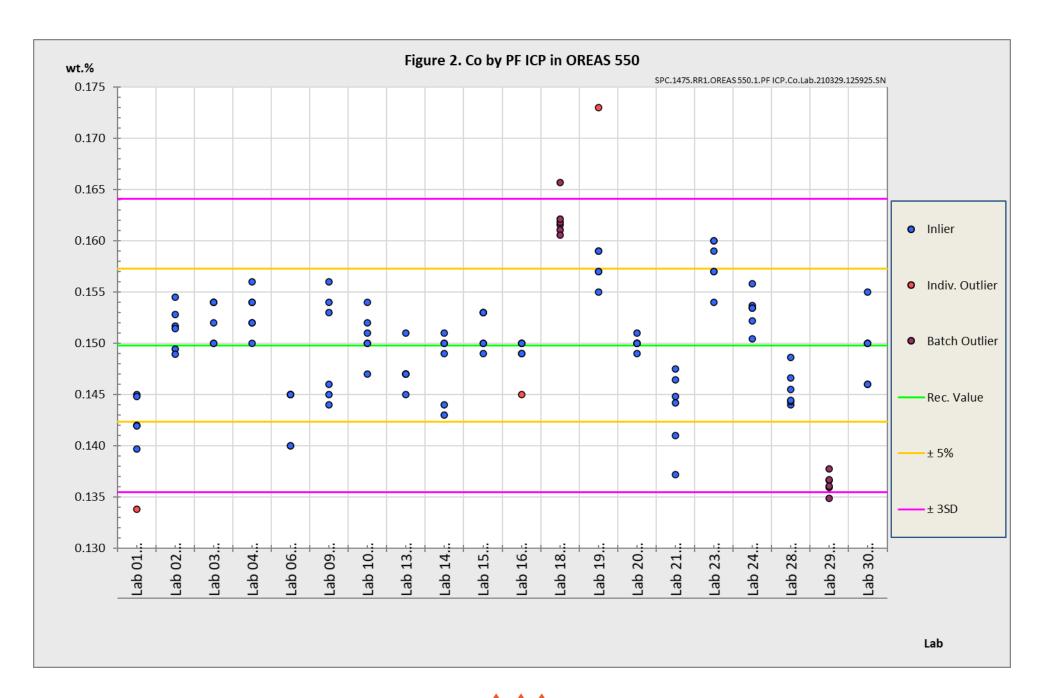

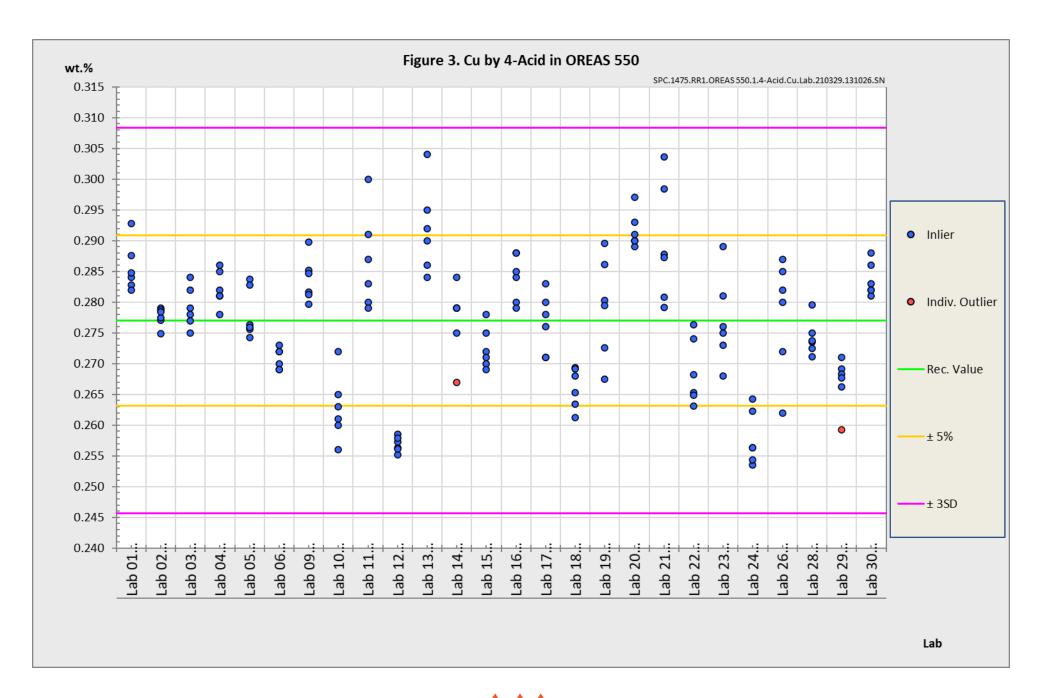

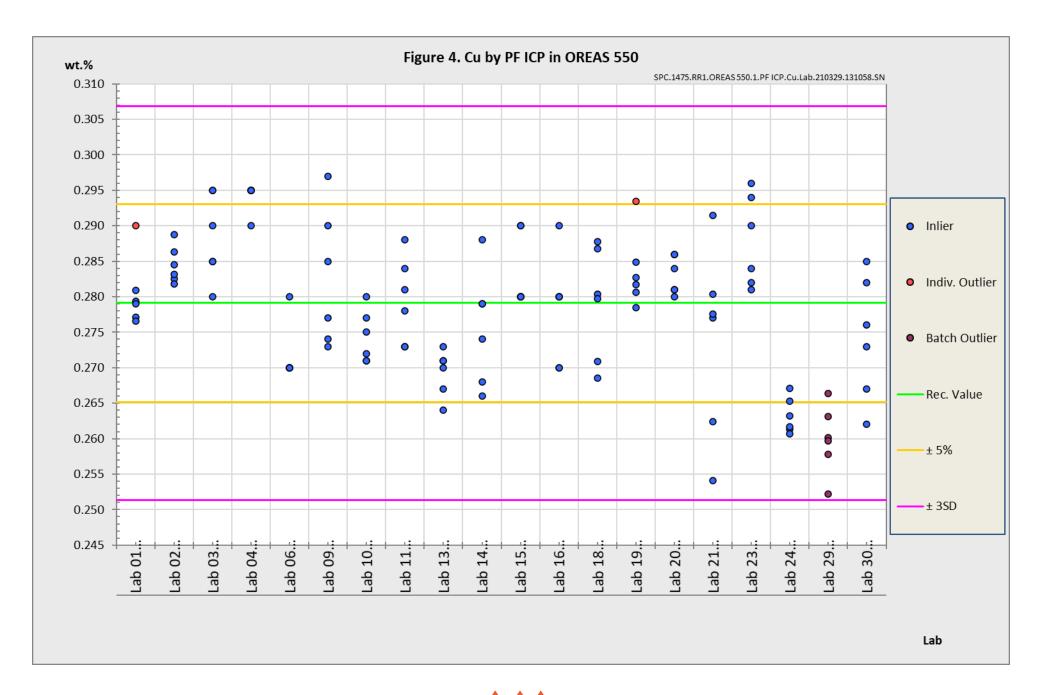

METROLOGICAL TRACEABILITY

The analytical samples were selected in a manner to represent the entire batch of prepared CRM. This 'representivity' was maintained in each submitted laboratory sample batch and ensures the user that the data is traceable from sample selection through to the analytical results that underlie the consensus values. Each analytical data set has been validated by its assayer through the inclusion of internal reference materials and QC checks during analysis.

The laboratories were chosen on the basis of their competence (from past performance in inter-laboratory programs undertaken by ORE Pty Ltd) for a particular analytical method, analyte or analyte suite, and sample matrix. Most of these laboratories have and maintain ISO 17025 accreditation. The certified values presented in this report are calculated from the means of accepted data following robust statistical treatment as detailed in this report.

Guide ISO/TR 16476:2016, section 5.3.1 describes metrological traceability in reference materials as it pertains to the transformation of the measurand. In this section it states, "Although the determination of the property value itself can be made traceable to appropriate units through, for example, calibration of the measurement equipment used, steps like the transformation of the sample from one physical (chemical) state to another cannot. Such transformations may only be compared with a reference (when available), or among themselves. For some transformations, reference methods have been defined and may be used in certification projects to evaluate the uncertainty associated with such a transformation. In other cases, only a comparison among different laboratories using the same method is possible. In this case, certification takes place on the basis of agreement among independent measurement results (see ISO Guide 35:2006, Clause 10)."

COMMUTABILITY

The measurements of the results that underlie the certified values contained in this report were undertaken by methods involving pre-treatment (digestion/fusion) of the sample. This served to reduce the sample to a simple and well understood form permitting calibration using simple solutions of the CRM. Due to these methods being well understood and highly effective, commutability is not an issue for this CRM. Being matrix-matched, OREAS 550 will display similar behaviour in the relevant measurement process to the routine 'field' samples for which OREAS 550 is designated to monitor. To maintain commutability, care should be taken to always ensure 'matrix matching' as close as practically achievable. The matrix and mineralisation style of the CRM is described in the 'Source Material' section and users should select appropriate CRMs matching these attributes to their field samples.

INTENDED USE

OREAS 550 is intended to cover all activities needed to produce a measurement result. This includes extraction, possible separation steps and the actual measurement process (the signal producing step). OREAS 550 may be used to calibrate the entire procedure by producing a pure substance CRM transformed into a calibration solution. OREAS 550 is intended for the following uses:

• For the monitoring of laboratory performance in the analysis of analytes reported in Table 1 in geological samples;

- For the verification of analytical methods for analytes reported in Table 1;
- For the calibration of instruments used in the determination of the concentration of analytes reported in Table 1.

STABILITY AND STORAGE INSTRUCTIONS

OREAS 550 has been prepared from a blend of copper-cobalt oxide ore samples. It is low in reactive sulphide (0.23 wt.% S) and in its unopened state and under normal conditions of storage has a shelf life beyond ten years. Its stability will be monitored at regular intervals and purchasers notified if any changes are observed.

Single-use sachets

Following analysis of the CRM subsample it is the manufacturers' expectation that any remaining material is discarded. The stability of the material after opening the sachet is not within the scope of proper use. However, if opened sachets are resealed after opening, then under ordinary* storage conditions the CRM will have a shelf-life beyond ten years.

*ordinary storage conditions: means storage not in direct sunlight in a dry, clean, well ventilated area at temperatures between -5° and 50°C.

Repeat-use packaging (e.g., 1kg plastic jars)

The stability of the CRM after opening the lid of the plastic jar is only affected by local atmospheric conditions with regard to oxidation and hygroscopic change. There is no segregation affect (please see our <u>Technical Note on Particle Segregation</u>).

The primary cause of change through oxidation is in relation to the breakdown of sulphide minerals to sulphates and is negligible for OREAS 550 given its low sulphur concentration (0.23 wt.% S).

Hygroscopic change is the amount of absorbed moisture (weakly held H₂O- molecules on the surface of exposed material) following exposure to the local atmosphere. Usually, equilibration of material to the local atmosphere will only occur if the material is spread into a thin (~2mm thick) layer and left exposed for a period of 2 hours. OREAS 550 contains a non-hygroscopic matrix and therefore, exposure to a local atmosphere that is significantly different (in terms of temperature and humidity) from the climate during manufacturing will have negligible impact on the precision of results. The 'Physical Properties' section indicates the approximate moisture concentration.

INSTRUCTIONS FOR CORRECT USE

The certified values for OREAS 550 refer to the concentration level in its packaged state. It should not be dried prior to weighing and analysis. 1kg jars permit repeated sampling as long as the lid is promptly re-secured to prevent airborne contamination.

Minimum sample size

As a practical guide, the minimum mass of sample used should match the typical mass that the laboratories used in the interlaboratory (round robin) certification program. This

means that different sample masses should be used depending on the operationally defined methodology.

- Lithium borate fusion with X-ray fluorescence finish: ≥0.2g;
- Loss on Ignition (LOI) at 1000°C: ≥1g;
- C and S by infrared combustion furnace/CS analyser: ≥0.1g;
- Cu and Co by 5% sulphuric acid leach with ICP or AAS finish: 0.5g;
- Cu and Co by 10% sulphuric acid leach with ICP or AAS finish: ≥0.5g;
- Sodium peroxide fusion with ICP-OES and/or MS finish: ≥0.2g;
- 4-acid digestion with ICP-OES and/or MS finish: ≥0.25g;
- Aqua regia digestion with ICP-OES and/or MS finish: ≥0.5g.

QC monitoring using multiples of the Standard Deviation (SD)

In the application of SD's in monitoring performance it is important to note that not all laboratories function at the same level of proficiency and that different methods in use at a particular laboratory have differing levels of precision. Each laboratory has its own inherent SD (for a specific concentration level and analyte-method pair) based on the analytical process and this SD is not directly related to the round robin program.

The majority of data generated in the round robin program was produced by a selection of world class laboratories. The SD's thus generated are more constrained than those that would be produced across a randomly selected group of laboratories. To produce more generally achievable SD's the 'pooled' SD's provided in this report include inter-laboratory bias. This 'one size fits all' approach may require revision at the discretion of the QC manager concerned following careful scrutiny of QC control charts.

The performance gates shown in Table 1 are intended only to be used as a first principle guide as to what a laboratory may be able to achieve. Over a period of time monitoring your own laboratory's data for this CRM, SD's should be calculated directly from your own laboratory's process. This will enable you to establish more specific performance gates that are fit for purpose for your application as well as the ability to monitor bias. If your long-term trend analysis shows an average value that is within the 95% confidence interval then generally there is no cause for concern in regard to bias.

It is important to note that in the analytical industry there is no standardisation of the aqua regia digestion process nor the 5% and 10% sulphuric acid leach process. These methods are partial empirical digests and differences in recoveries for various analytes are commonplace. These are caused by variations in the digest conditions and in the case for aqua regia digestion can include the ratio of nitric to hydrochloric acids, acid strength, temperatures, leach times and secondary digestions. Recoveries for sulphide-hosted base metal sulphides approach total values, however, other analytes, in particular the lithophile elements, show greater sensitivity to method parameters. This can result in lack of consensus in an inter-laboratory certification program for these elements.

The approach applied here is to report certified values in those instances where reasonable agreement exists amongst a majority of participating laboratories. The results of specific laboratories may differ significantly from the certified values, but will, nonetheless, be valid and reproducible in the context of the specifics of the aqua regia method in use. Users of this reference material should, therefore, be mindful of this limitation when applying the certified values in a quality control program.

HANDLING INSTRUCTIONS

Fine powders pose a risk to eyes and lungs and therefore standard precautions such as the use of safety glasses and dust masks are advised.

LEGAL NOTICE

Ore Research & Exploration Pty Ltd has prepared and statistically evaluated the property values of this reference material to the best of its ability. The Purchaser by receipt hereof releases and indemnifies Ore Research & Exploration Pty Ltd from and against all liability and costs arising from the use of this material and information.

DOCUMENT HISTORY

Revision No.	Date	Changes applied
1	21 st April, 2021	Applied minor statistical correction to Cobalt by XRF method results.
0	7 th April, 2021	First publication.

QMS ACCREDITATION

ORE Pty Ltd is accredited to ISO 9001:2015 by Lloyd's Register Quality Assurance Ltd for its quality management system including development, manufacturing, certification and supply of CRMs.





CERTIFYING OFFICER

Sp

21st April, 2021

Craig Hamlyn (B.Sc. Hons - Geology), Technical Manager - ORE P/L

REFERENCES

Govett, G.J.S. (1983), ed. Handbook of Exploration Geochemistry, Volume 2: Statistics and Data Analysis in Geochemical Prospecting (Variations of accuracy and precision), P.O. Box 330, 1000 AH Amsterdam, The Netherlands.

ISO Guide 30:2015, Terms and definitions used in connection with reference materials.

ISO Guide 31:2015, Reference materials – Contents of certificates and labels.



ISO Guide 35:2017, Certification of reference materials - General and statistical principals.

ISO 16269:2014, Statistical interpretation of data – Part 6: Determination of statistical tolerance intervals.

ISO/TR 16476:2016, Reference Materials – Establishing and expressing metrological traceability of quantity values assigned to reference materials.

ISO 17025:2005, General requirements for the competence of testing and calibration laboratories.

Munsell Rock Color Book (2014). Rock-Color Chart Committee, Geological Society of America (GSA), Minnesota (USA).

APPENDIX

For Cu and Co by the two sulphuric acid leaches, specific methodologies were detailed for the participating laboratories to follow:

5% sulphuric acid leach

- 1. Weigh 0.500 ± 0.002g of sample pulp into a clean 250 ml flask.
- 2. Add to the flask 0.5g of Sodium Sulphite (AR Grade).
- 3. Add 50 ml of approximately 55 g/l Sulphuric acid solution (prepared from a 98% pure concentrated sulphuric acid).
- 4. Put the cap on the flask and start automatic shaking. Leave the sample on continuous shaking to leach for four (4) hours.
- 5. Remove the cap and add 25 ml of hydrochloric acid.
- 6. Dilute to a final volume of 250 ml with distilled/deionised water up to the mark and mix again by inverting at least 10 times. Allow the solution to settle for 30 minutes.
- 7. The solution is now ready to be analysed by ICP or AAS.

10% sulphuric acid leach

- 1. Weigh 0.500 ± 0.002g of sample pulp into a clean 250 ml flask.
- 2. Add to the flask 0.5g of Sodium Sulphite (AR Grade).
- 3. Add 50 ml of approximately 110 g/l Sulphuric acid solution (prepared from a 98% pure concentrated sulphuric acid).
- 4. Put the cap on the flask and start automatic shaking. Leave the sample on continuous shaking to leach for four (4) hours.
- 5. Remove the cap and add 25 ml of hydrochloric acid.
- 6. Dilute to a final volume of 250 ml with distilled/deionised water up to the mark and mix again by inverting at least 10 times. Allow the solution to settle for 30 minutes.
- 7. The solution is now ready to be analysed by ICP or AAS.