


DigiPREP CUBE

The **DigiPREP CUBE** digestion system is an innovative, 25-position digestion block employing 16 mm tubes, designed for Chemical Oxygen Demand and other digestions. Simple operations can be achieved through keypad selection; by setting the required time and temperature, or by selecting a predefined program.



- Offers 3-predefined programs
 - COD DIN @ 148°C in 120 minutes
 - Fast COD @ 150°C in 15 minutes
 - COD EPA @ 150°C in 120 minutes
- Stores 2-user defined programs for added flexibility
 - Temperature: Ambient to 180°C
 - Time: 0 to 360 minutes
- Useable in 7 different languages

Special Promotion !!
 With the purchase of a
DigiPREP CUBE receive a **FREE**
 box of **COD Digestion Tubes**.

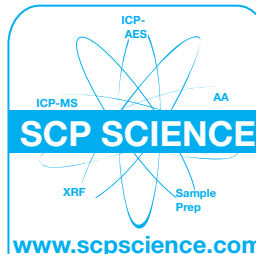

Sample Prep
ICP-AES & MS
AA
XRF
Calibration Standards
Certified Reference Materials

Ordering Information & Specifications

DigiPREP CUBE


Catalog Number (115V)	010-510-003	Over-Temp Protection	Yes
Catalog Number (230V)	010-510-004	Heating Block	Teflon-coated graphite
Capacity	25 Tubes	Power (010-510-003)	115V / 3A
Temperature	Ambient to 180°C	Power (010-510-004)	230 / 1.5A
Timer	0 to 360 Minutes	Size (Inches)	8W x 8D x 8.75H
	Digital display to 1 min accuracy	(cm)	20.5W x 20.5D x 23H
	Audible end-of-cycle alarm	Certification	CE/UL/CSA
Set-Temperature	Precision within 0.1°C of Setpoint		

Options and Accessories

Description	Concentration	Quantity	Catalog #
COD Digestion Tubes	100 ppm O ₂	25 + 1 QC	250-130-006
COD Digestion Tubes	1000 ppm O ₂	25 + 1 QC	250-130-016
COD Digestion Tubes	10 000 ppm O ₂	25 + 1 QC	250-130-026

Description	Concentration	20 ml	125 ml	500 ml	1 L
COD Standard	100 ppm O ₂	---	250-130-512	250-130-550	250-130-551
COD Standard	1000 ppm O ₂	---	250-130-602	250-130-600	250-130-601
COD Standard	10 000 ppm O ₂	---	250-130-652	250-130-650	250-130-651
CRM for COD	Low Level (200 ppm)	140-715-101	---	---	---
CRM for COD	High Level (1000 ppm)	140-715-102	---	---	---

COD Digestion Tubes

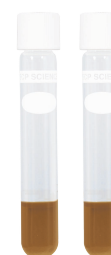

AccuSPEC COD Digestion Tubes use standard 16 mm tubes that have been prefilled with digestion reagent.

- Available in 3 ranges
- Compatible with standard reactors and spectrometers
- Fully compatible with EPA 410.4 and Standard

MK-MKG003-DGCB-1.0-E

For more information concerning the **DigiPREP CUBE** or options & accessories, please contact your sales representative or distributor.

USA
 Tel.: (800) 361-6820
 Fax: (800) 253-5549

Canada / International
 Tel.: (800) 361-6820 / (514) 457-0701
 Fax: (800) 253-5549 / (514) 457-4499

Europe
 Tel.: +33 (0)1 69 18 71 17
 Fax: +33 (0)1 60 92 05 67